




Ksure®

Silicone Foam Dressing

A trusted solution for effective wound management and healing

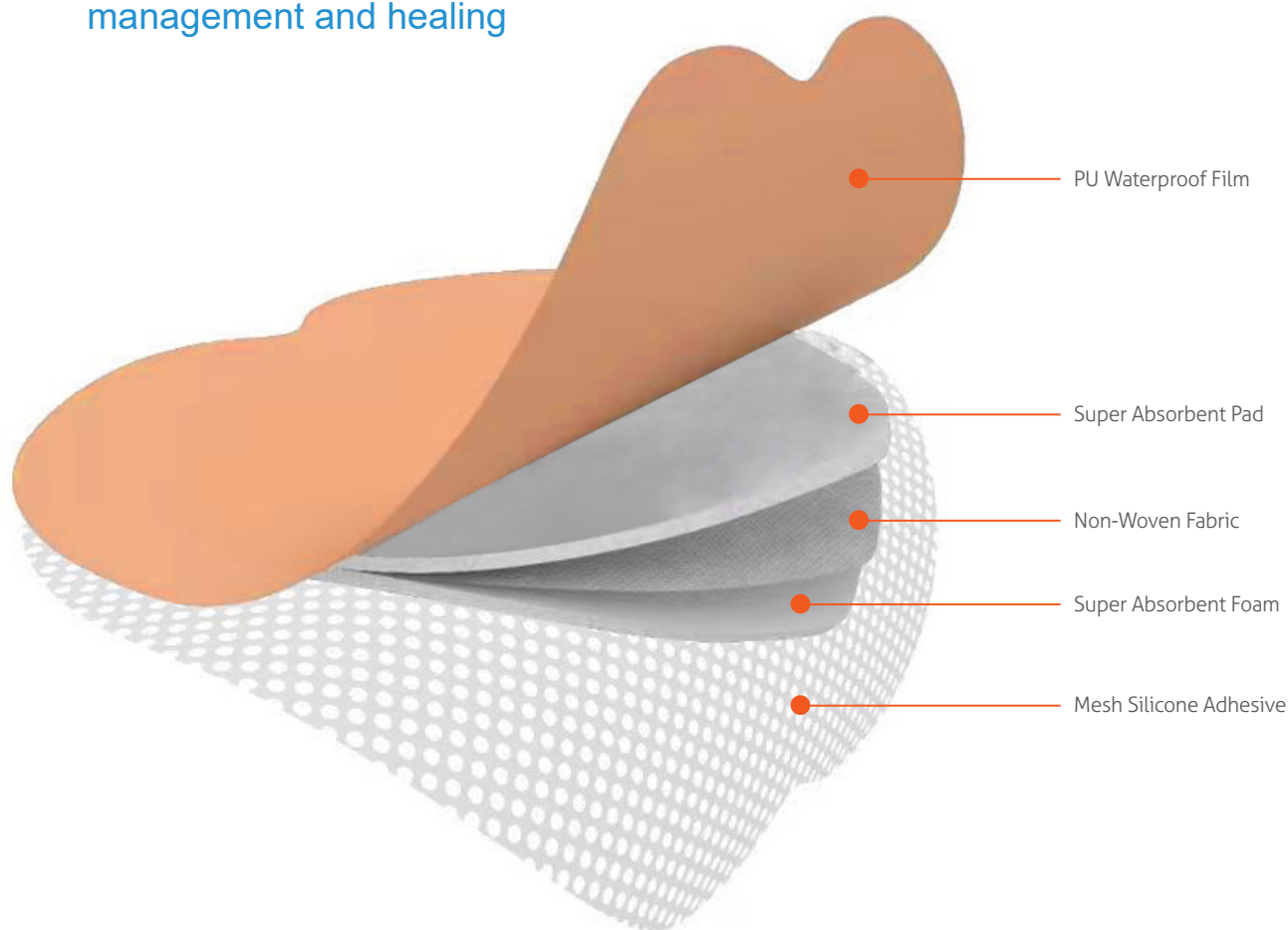

KB MEDICAL(GROUP),INC.

Silicone Foam Dressing

Silicone Foam Dressing combines the advantages of polyurethane foam, such as softness, skin-friendliness, elasticity and vertical absorption of exudate, with the advantages of silicone, including low tackiness, safety, non-irritation, and reusability. They exhibit excellent performance in terms of exudate absorption, absorption rate, application viscosity, emulsion ratio, elasticity coefficient, and softness.

Bordered Silicone Foam Dressing

Bordered Silicone Foam Dressing is with a five-layer structure:

1. PU Waterproof Film
2. Super Absorbent Pad
3. Non-Woven Fabric
4. Super Absorbent Foam
5. Mesh Silicone Adhesive



Feature

- Strong absorbency: Suitable for wounds with hypertonic exudate
- Bordered design: Provides additional fixation force
- Gentle Silicone Adhesive: Ideal for delicate skin.
- Bacterial and Viral Barrier
- Long Wear Time: Can often remain in place for up to 7 days

Application

The Bordered Silicone Foam Dressing is designed for a wide range of exuding wounds such as pressure ulcers, leg and foot ulcers, traumatic wounds (e.g. skin tears) and surgical wounds. Bordered Silicone Foam Dressing reduces postoperative blistering and may also be used as part of a prophylactic therapy to help prevent skin damage, e.g. pressure ulcers.


Ordering information

Product Size		Pcs/Box	Box/Case	Code	Product Size		Pcs/Box	Box/Case	Code
Inch	mm				Inch	mm			
3 x 3	75 x 75	10	20	SFB503030	6 x 8	150 x 200	10	10	SFB506080
4 x 4	100 x 100	10	20	SFB504040	8 x 8	200 x 200	5	10	SFB508080
4 x 6	100 x 150	10	10	SFB504060	5 x 8	125 x 200	10	10	SFB505080
4 x 8	100 x 200	10	10	SFB504080	7.3 x 9.5	185 x 240	10	5	SFB507395
4 x 10	100 x 250	10	10	SFB5040100	6 x 7.1	150 x 180	10	20	SFB506071
4 x 12	100 x 300	10	10	SFB5040120	7.1 x 7.1	180 x 180	5	10	SFB507171
5 x 5	125 x 125	10	20	SFB505050	9 x 9	225 x 225	5	10	SFB509090
6 x 6	150 x 150	10	10	SFB506060	*Can be customized				

Silicone Foam Dressing Without Border

Silicone Foam Dressing Without Border is with a three-layer structure:

1. PU Waterproof Film
2. Super Absorbent Foam
3. Mesh Silicone Adhesive

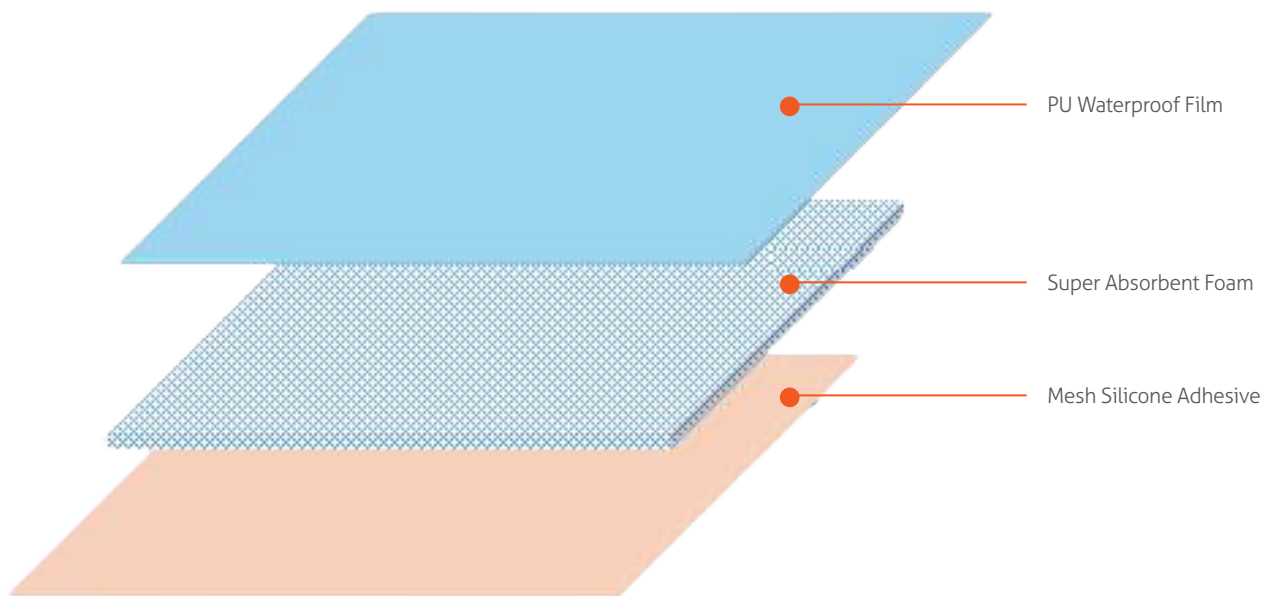

Feature

- Reduces wound pain.
- Minimizes damage to surrounding skin.
- Prevents wound infection and decay.
- Keeps the wound moist, which helps accelerate wound healing.

Application

The Silicone Foam Dressing Without Border is suitable for secondary healing of superficial granulating wounds, chronic and acute exudative wounds, and full-thickness or partial-thickness wounds (such as pressure ulcers, leg ulcers, diabetic foot ulcers, infected wounds, malignant lesions, first and second-degree burns, donor sites, skin tears, fungal ulcers, and other types of wounds). Wounds with necrosis or eschar can also be treated in combination with gel dressings.

Silicone Foam Dressing Without Border, 3 Layer 5mm

5mm thickness design enhances the absorption capacity and reduces pressure through buffering.

Ordering information

Product Size		Thickness (mm)	Pcs/Box	Box/Case	Code
Inch	mm				
2 x 2	50 x 50	5	10	20	SF352020
3 x 3	75 x 75	5	10	20	SF353030
4 x 4	100 x 100	5	10	20	SF354040
4 x 6	100 x 150	5	10	10	SF354060
4 x 7	100 x 175	5	10	10	SF354070
4 x 8	100 x 200	5	10	20	SF354080
4 x 10	100 x 250	5	10	10	SF3540100
4 x 12	100 x 300	5	10	10	SF3540120
5 x 5	125 x 125	5	10	20	SF355050
6 x 6	150 x 150	5	10	10	SF356060
5.5 x 8.3	140 x 210	5	10	10	SF355583
6 x 7.1	150 x 180	5	10	10	SF356071

Silicone Foam Dressing Without Border, 3 Layer 2mm

2mm thickness design is suitable for joints and frequently moving body parts.

Ordering information

Product Size		Thickness (mm)	Pcs/Box	Box/Case	Code
Inch	mm				
2 x 2	50 x 50	2	10	20	SF322020
3 x 3	75 x 75	2	10	20	SF323030
4 x 4	100 x 100	2	10	20	SF324040
4 x 6	100 x 150	2	10	20	SF324060
4 x 7	100 x 175	2	10	10	SF324070
4 x 8	100 x 200	2	10	10	SF324080
4 x 10	100 x 250	2	10	20	SF3240100
4 x 12	100 x 300	2	10	10	SF3240120
5 x 5	125 x 125	2	10	20	SF325050
6 x 6	150 x 150	2	5	10	SF326060
5.5 x 8.3	140 x 210	2	10	10	SF325583
6 x 7.1	150 x 180	2	10	10	SF326071

Post-Op Silicone Dressing

Post-Op Silicone Dressing is advanced wound care products designed for surgical wounds to manage exudate, reduce pain during dressing changes, and promote a moist healing environment.

Post-Op Silicone Dressing is with a three-layer structure:

1. PU Waterproof Film
 2. Super Absorbent Pad
 3. Mesh Silicone Adhesive
- * Can be customized.



Feature

- Specifically designed for post-operative recovery.
- Strong absorbency: Improved exudate management and enhanced absorption efficiency.
- Bordered design: Provides additional fixation force.
- Gentle Silicone Adhesive: The silicone layer is gentle during removal, reducing secondary trauma.
- Infection Barrier: Protects the surgical site from external contaminants.

Applications

The Post-Op Silicone Dressing is mainly intended for managing acute surgical wounds by absorbing exudate while providing a non-adherent, protective barrier. They are specifically engineered to maintain a moist environment that accelerates healing and reduces the risk of scarring.

Ordering information

Product Size		Pcs/Box	Box/Case	Code
Inch	mm			
3 x 3	75 x 75	10	20	SPO303030
4 x 4	100 x 100	10	20	SPO30440
4 x 6	100 x 150	10	10	SPO304060
4 x 8	100 x 200	10	10	SPO304080
4 x 10	100 x 250	10	10	SPO3040100
4 x 12	100 x 300	10	10	SPO3040120
5 x 5	125 x 125	10	20	SPO305050
6 x 6	150 x 150	10	10	SPO306060
6 x 8	150 x 200	10	10	SPO306080
8 x 8	200 x 200	5	10	SPO308080
5 x 8	125 x 200	10	10	SPO305080
7.3 x 9.5	185 x 240	10	5	SPO307395
6 x 7.1	150 x 180	10	20	SPO306071
7.1 x 7.1	180 x 180	5	10	SPO307171
9 x 9	225 x 225	5	10	SPO309090

Flex Silicone Foam Dressing

Flex Silicone Foam Dressing is an advanced, multi-layered wound dressings featuring proprietary "Flex Technology" for enhanced adaptability and secure fit, especially on mobile areas or difficult body contours.

Flex Silicone Foam Dressing is with a five-layer structure:

1. PU Waterproof Film
2. Super Absorbent Pad
3. Non-Woven Fabric
4. Super Absorbent Foam
5. Mesh Silicone Adhesive



Feature

- High Flexibility: Pre-cut design in the dressing's spreading and retention layers, allowing for a 360-degree stretch
- Enhanced Conformability: Can stretch in all directions helps it conform well to challenging anatomical areas like joints (heels, elbows, knees), ensuring it stays securely in place for longer periods.
- Strong absorbency: Improved exudate management and enhanced absorption efficiency.
- Bordered design: Provides additional fixation force.
- Gentle Silicone Adhesive: The silicone layer is gentle during removal.
- Infection Barrier: Protect the surgical site from external contaminants.
- Long Wear Time: Can often remain in place for up to 7 day.

Ordering information

Product Size		Pcs/Box	Box/Case	Code
Inch	mm			
3 x 3	75 x 75	10	20	SFBF503030
4 x 4	100 x 100	10	20	SFBF504040
4 x 6	100 x 150	10	10	SFBF504060
4 x 8	100 x 200	10	10	SFBF504080
4 x 10	100 x 250	10	10	SFBF5040100
4 x 12	100 x 300	10	10	SFBF5040120
5 x 5	125 x 125	10	20	SFBF505050
6 x 6	150 x 150	10	10	SFBF506060
6 x 8	150 x 200	10	10	SFBF506080
8 x 8	200 x 200	5	10	SFBF508080
5 x 8	125 x 200	10	10	SFBF505080
7.3 x 9.5	185 x 240	10	5	SFBF507395
6 x 7.1	150 x 180	10	20	SFBF506071
7.1 x 7.1	180 x 180	5	10	SFBF507171
9 x 9	225 x 225	5	10	SFBF509090

Applications

The Flex technology provides a robust, yet comfortable, solution for effective wound care, especially for active patients or wounds in difficult-to-dress locations.



Scan QR code to learn more.

KB MEDICAL (GROUP), INC. Ksure® is a registered brand of KB Medical

Add.: 5325 Alton Parkway Ste C 732, Irvine, CA 92604 Tel: 1-949-333-3785 Email: info@kbmedcare.com

Website: www.kbmedcare.com ; www.kbmed.com

